



GOOD DESIGN AWARD 2015
GRAND AWARD



reddot award 2014
honourable mention

IDEX
INTERNATIONAL DESIGN
EXCELLENCE AWARDS '14
BRONZE





WHILL

Design your own road

It began with one user's voice.

"I've even given up on going to the grocery store two blocks away."

Our journey started with this simple statement from a wheelchair user. He did not want to be seen in public as a wheelchair user. He disliked the negative perceptions associated with wheelchairs; that he must be ill or weak.

To solve this, we started to create an original mobility device just for him. We exhibited our first prototype at the Tokyo Motor Show, where we were overwhelmed by the positive reaction we received from all over the world.

We realized that many people desired this kind of innovation in the mobility industry.

We saw how meaningful what we were creating could be for people with mobility issues.

Nobody should be confined because they don't feel comfortable with how they move from point A to point B.

We want every individual to be confident, comfortable, and independent.

This is why we founded WHILL.



Meet WHILL

WHILL Model M doesn't hold you back, it powers you forward.

With simplistic beauty and revolutionary design, Model M provides mobility in a versatile, stylish and unobtrusive way.



4-Wheel Drive

Live your life without disruption and go wherever you want to go.



Long-Distance Range

Worry-free extended range of up to 15.12 miles to accomplish all your daily tasks.



Obstacle Clearance

Large 9.8" front wheels will clear terrain obstacles in your way.



All Directional Wheels

Drive the highest performance powered wheelchair in the industry.



Climbing Capability

Climb steep hills up to 10° incline and take on any environment.



Additional Support

Model M comes with various back support options and many more optional additions to fit your body type.



Design

High performance meets unparalleled beauty

WHILL is an extension of a life beautifully lived. Its unparalleled navigational and functional capabilities, ease of use, compact form and clean design provide freedom, independence and individualized expression that puts your life front and center.



Compact 23.6 inch wide frame



Simple user interface



Durable support system



Moveable arms for easier transfers


Seat slides up to 5.9 inches forward



Technology


Maneuver the world like never before

WHILL was founded and designed by former Toyota and Nissan automobile industry experts and engineers. It is the world's first powered wheelchair to feature All Directional Wheel Technology combined with 4-wheel drive, a proprietary core technology that enables advanced motion and multi-terrain coverage. That means WHILL can take you from your home to your urban office, down the street to the cafe and even up to the mountains for vacation.

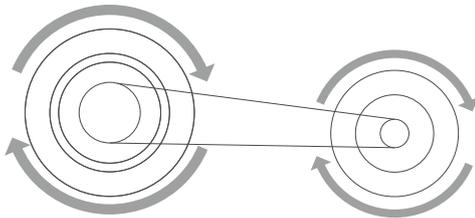

4WD

Both front and rear wheels move together to provide maximum traction on any terrain.

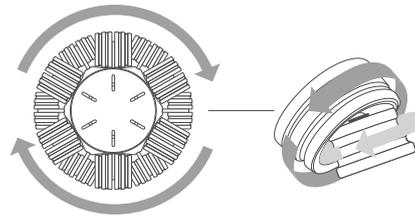


All Directional Wheel

24 rollers on each wheel allow WHILL to turn on its own axis



WHILL turns on its own axis with a 28 inch turning radius

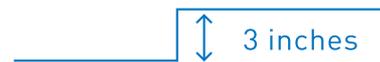

Clear obstacles of up to 3 inches.



Climb inclines of up to 10 degrees



Easily drive over rough and loose surfaces, including gravel.



Configurations

Mobility that meets your standards

Model M can be customized to fit your body type and support needs. Do you want the Varilite Icon Back System in Short to keep a low profile? Or would you prefer the Varilite Icon Back System in Deep for a more encompassing fit? Either way, WHILL has an option for you.



Varilite Icon Back System
Short



Varilite Icon Back System
Mid



Varilite Icon Back System
Tall



Varilite Icon Back System
Deep



Body Point U-Shaped
Joystick Handle



Body Point EvoFlex Lap
Belt
(Included)



Pressure Relief Handles
(Included)



Varilite Meridian Seat
Cushion
(Included)



Rear Tie Down Points
(Included)

Accessorize

Have it all, or not. It only depends on what you want.


Build WHILL your way by choosing only the accessories that matter to you. These are some of the options that we can ship with your device. Many of them are also available from third party sellers.



WHILL comes in
either white or black



Skip Hop Organizer

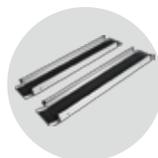

WheelzAhead Portable
Ramps



KNOG Headlights




WHILL Cover


Fitting

Fully configurable seating and positioning allows for a comfortable and personalized fit.

Everything from the back support height to the foot support length is adjustable. We fit it to your measurements before shipping, and someone is always there to make sure its perfect on the day of delivery.

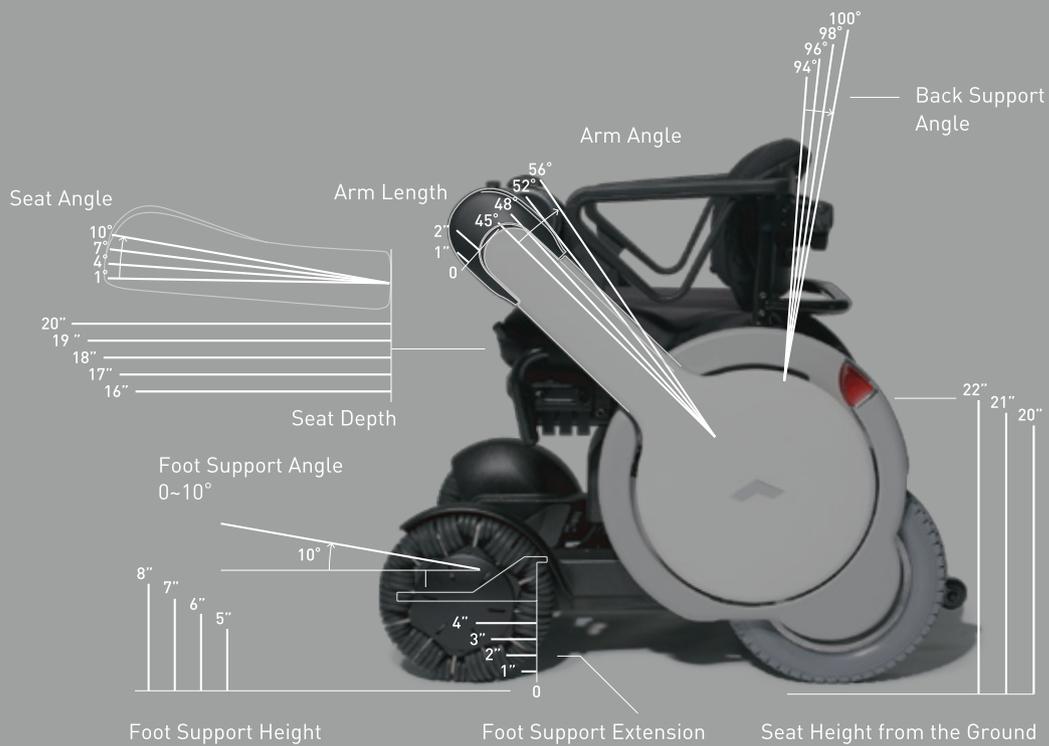


The directional controller can be on the right or left

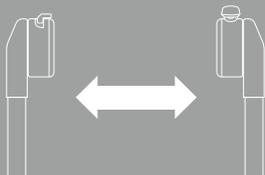

Back Support

The back of Model M is a simple pipe structure with two brackets set on each side. The brackets are compatible with the Varilite VariLock Mounting Kit which are compatible with the Icon Back System.

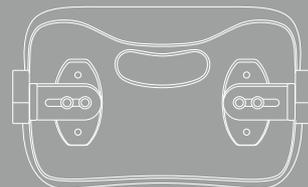

How to Purchase

Contact Us


1. Fill out the Test Drive Request form on our website, or call (844) 699-4455 extension 2 to request a test drive.
2. We will help you locate your nearest reseller, who will then schedule the test drive at your earliest convenience.

Meet with your physician

3. Talk to your physician about prescribing WHILL for your personal mobility needs. Your doctor will determine the medical necessity of providing you with a powered wheelchair and seating evaluation, as evidenced by the prescription.
4. Once payment is received, your order can be placed. Your reseller will fill out our order form with your information including body measurements, so we can begin to fit the chair to you before it even ships.

Delivery and Training

5. A technician from your reseller location will be on site to deliver Model M to a location of your choice.
6. At the time of delivery, your technician will perfect the seating and train you on how to use your new powered wheelchair.

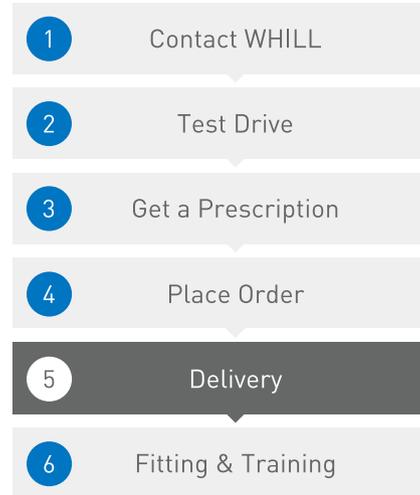

Over 40 Locations Nationwide

We've partnered with resellers from CA to NY with plenty of stops in between, and we're working quickly to bring WHILL to new locations each month.



WHILL Model M

WHILL moves seamlessly throughout your life and goes effortlessly wherever you want to go. It is the combination of both high performance and unparalleled beauty.

Spec Detail

Height	Depends on choice of back support
Width	23.6 in
Length	44.3" *
Weight	255 lb
Range	15.12 miles **
Power	2 x 12V 50Ah batteries
Speed	5.5 mph top speed
Climb	10° Incline Capability
Front Wheel	9.8" diameter
Back Wheel	12.5" diameter
Ground Clearance	3.5"
Weight Capacity	220 lb

*Figures shown in maximum setting

** Continuous Theoretical Driving Range





WHILL, Inc.

285 Old County Rd. #6 San Carlos, CA 94070

+1-844-MY-WHILL

whill.us

320-01452-MKT-01

03/16



WHILL Model M CAUTION:

Federal law restricts this device to sale by or on the order of a physician.